


**Product Name :**  
Soil Cone Penetromter

**Product Code :**  
LIM-CAT-L0044-00012

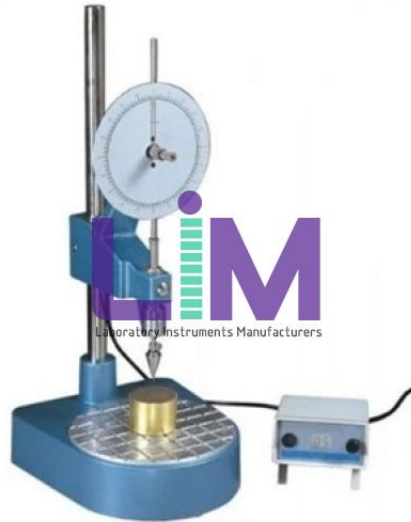

**Description :**

Soil Cone Penetromter

**Technical Specification :**

Soil Cone Penetromter

Soil Cone Penetrometer of Soil is based on the relationship between moisture content and the penetration of the cone into the soil sample. The moisture content at a given depth of penetration of the cone is taken as the liquid limit.

The Automatic Penetrometer (Digital Timer) is supplied complete with a stainless steel penetration test cone 35mm long with a smooth polished surface and an angle of 30° and a Penetration test cup 55mm dia x 40mm deep. Facility exists for adjusting cone height in relation to the specimen.

The apparatus is suitable for operation on 220V, Single Phase, 50Hz, AC Supply.

Stainless steel penetration test cone 35mm long.

Angle: 30°

Penetration test cup: 55mm x 40mm.

